

AdiaX

Prep

Device for the preparation of faeces samples

- The easiest sample preparation for ADIAVET™ PARATB REAL TIME
- One device
- Time saving (4X faster)



6 elements



1 element



All-in-one sample preparation device including:

- ✓ No weighing
- ✓ No pipetting
- ✓ No clogging

Easy to use for large amount of samples analysis in a context of herd monitoring

New procedure for easier sample preparation



Fastest sample preparation on the market:

- Only 1 spoon
- Buffer included (30 mL)



The most precise sample preparation on the market:

- High quantity of faeces (4,5 g)



The most practical sample preparation on the market:

- High filter preventing clogging
- Screw cap preventing splashing when closing
- Screw cap opens in 2 turns with one hand
- Double integral cap
- Standard diameter 50 mL tube: easy transport for 24 series, held on standard rack
- Tightness during transport



1



Sampling

2



Mix


3




Dispense

Easier sample prep : pre-processing of 96 samples in 1h with ADIAPREP™

60 min

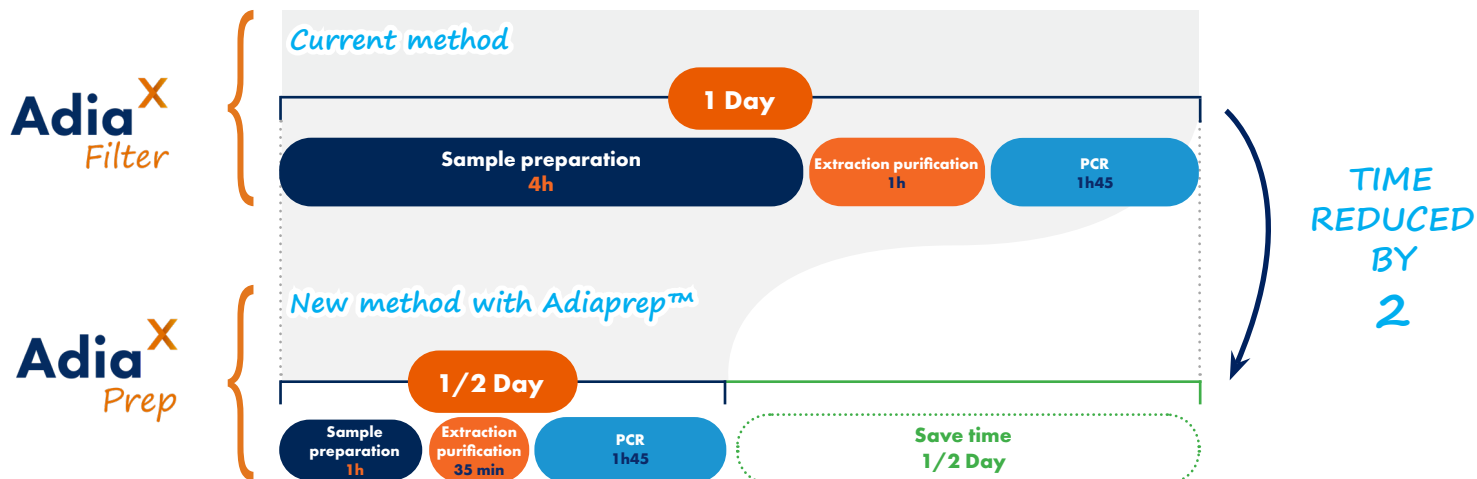

Sample preparation with ADIAPREP™

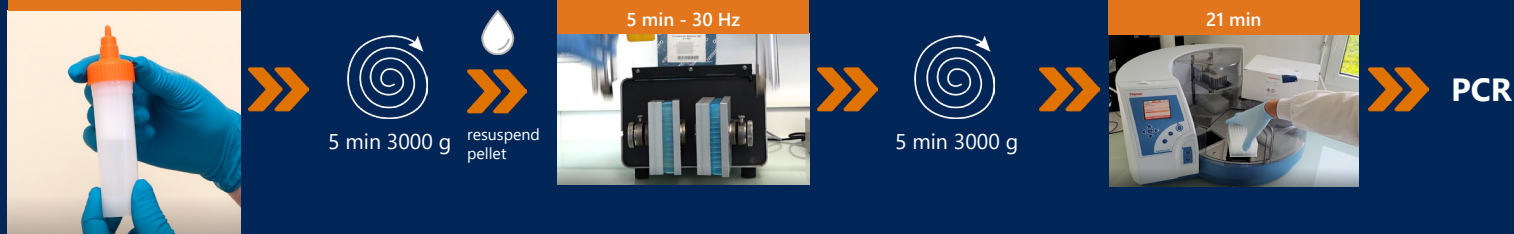

External validation with the LNCR

(EAVLD Sevilla communication Oct 2022 : Evaluation of new device and procedure for the detection and quantification of MAP Q_{PCR} on large series of samples)

- Easy to use with current material
- Large series of samples
- High correlation with internal methods
- Ready to use material



		Routine Qualified Samples	
		+	-
ADIAPREP™ + ADIAVET PARATB PCR	+	16	1
	-	1	23

TO PLACE AN ORDER :

Code	Description	Nb. of reactions
ADI045-100	ADIAVET™ PARATB REAL TIME	100 R
ADPBIAR-4x96	ADIAPURE™ GLASS BEADS RACKS 4x96	384 R
NADI003	ADIAMAG™ (magnetic beads extraction kit)	200 R
ADPREP-200	ADIAPREP™	200 R

ADDITIONAL KITS FOR METHOD ADOPTION AND qPCR :

Code	Description	Nb. of reactions
ADC04EPC	Extraction Positive Control PARATB	100 R
ADC04SQ01	Quantified Extraction Positive Control PARATB faeces	100 R
ADC04LD	LD PCR Positive Control - PARATB	100 R

Why ADIAPREP™ ?



Contact us

support.pcr@biox.com

+33 02 96 68 40 20

www.biox.com



Smart solutions for sharp decisions